



FITTING GUIDE

**Mounts & More Mini Floorstand Base
with Post Release**

MINF-SA

MINF-PA

MINF -FA8

Mounts & More Ltd

Fitting Guide

Mounts & More Mini Floorstand Base
with Post Release



Mounts & More Mini Floorstand Base with Post Release

Index:

<i>Page 2</i>	<i>Index</i>
	<i>1. stability advice</i>
<i>Page 3-4</i>	<i>2. The parts included in a MINF-SA, MINF-FA8 and MINF-PA</i>
<i>Page 5</i>	<i>3. Using the Post Release</i>
	<i>4. Using the Repositionable Height Adjust Handle</i>
<i>Page 6</i>	<i>5. Adjusting the Short Arm and balljoints</i>
	<i>6. Adjusting the Post Articulating Quick Release Base</i>

1. Stability

Although designed with safety in mind and remarkably stable for their size, these are small footprint bases and must be used responsibly.

- Do not exceed the prescribed weights (2kg for the MINF-SA, 4kg for the MINF-FA8 and 4kg for the MINF-PA).
- Keep the arms pointing straight over the centre of the base (see Fig. 1.1, 1.2 and 1.3).

Tip: Leaving the front castors unlocked where possible (especially appropriate on the MINF-PA) makes it even more unlikely that the Mini Floorstand can be accidentally tipped over.

Fig. 1.1



MINF-SA
Mini Floorstand with
Short Arm

Fig. 1.2




MINF-FA8
Mini Floorstand with
Fixed 8" Arm

Fig. 1.3




MINF-PA
Mini Floorstand with
Post Articulating Quick
Release Base

2. The parts included in the **Mounts & More Mini Floorstands** are:

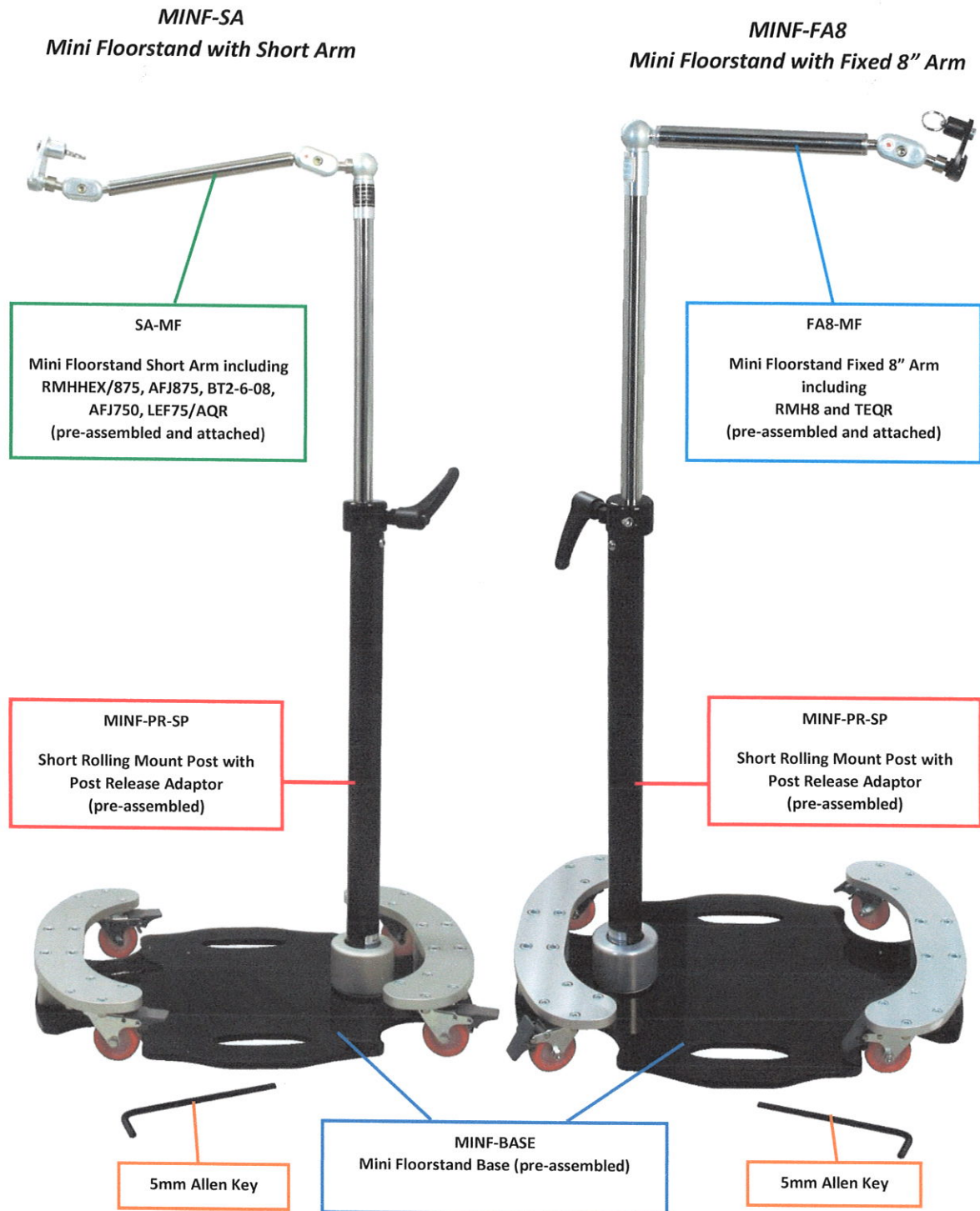


MINF-PA
Mini Floorstand with Post Articulating Quick
Release Base

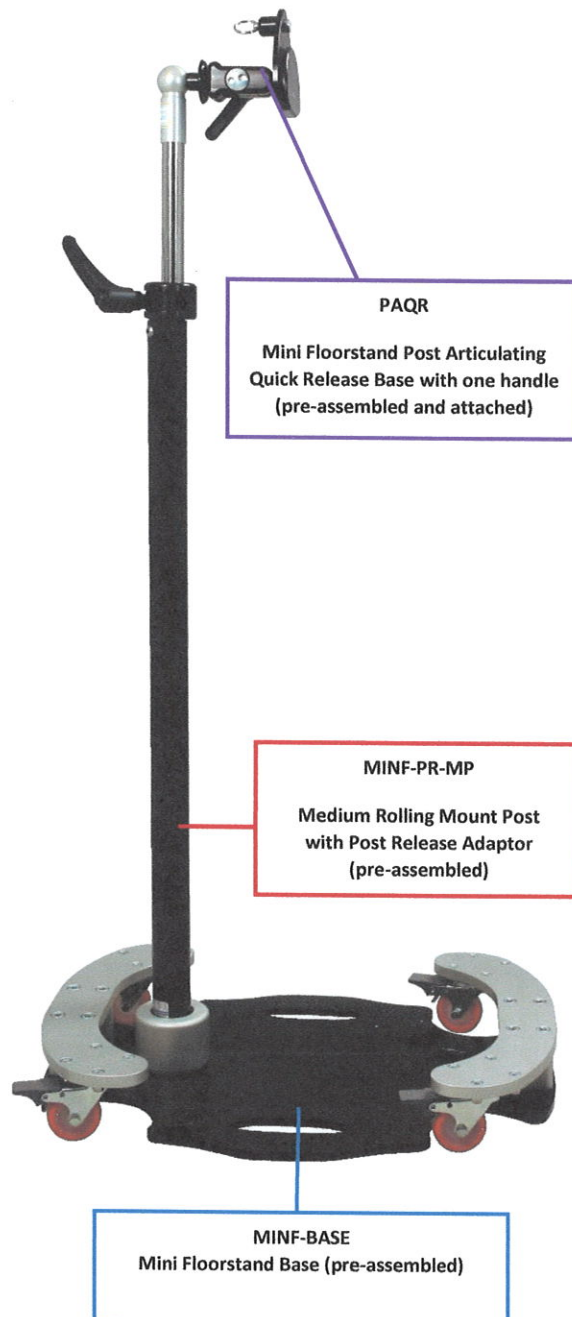


3. Attaching the **MINF-PR-SP** or **MINF-PR-MP** Post with Post Release Adaptor to the Base

The Post will only attach to the base one way. On the Base Receiver is a **peg** (see Fig. 3.1) which engages with a **hole** on underside of the Post Adaptor (see Fig. 3.2)

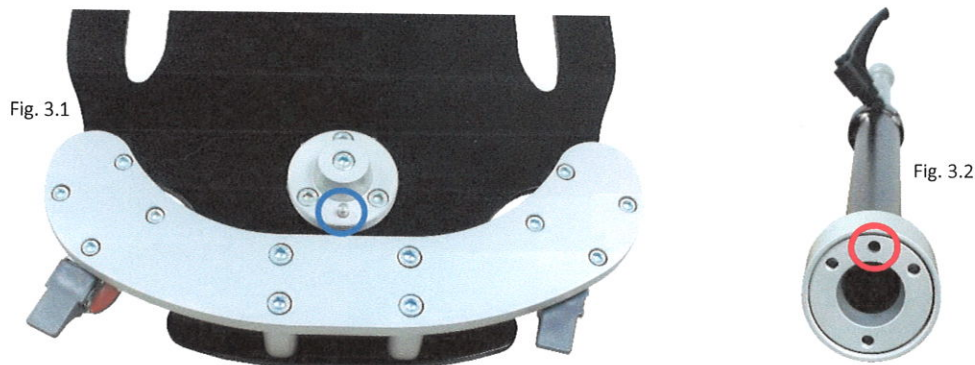

Loosen the height adjust handle (see Fig. 3.3) or raise the post tube then place the Post on the Base Receiver ensuring that the peg is properly engaged in hole in the Post Adaptor. Tighten the Post Release Collar (see Fig. 3.4).

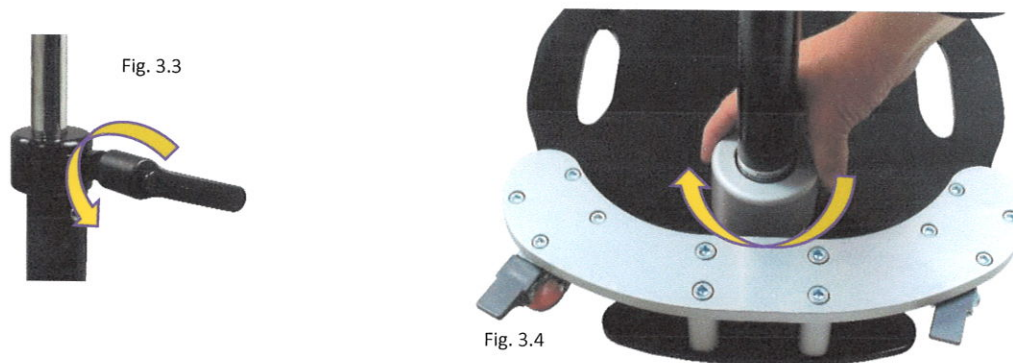


4. Adjusting the height of the Post tube

ALWAYS SUPPORT OR REMOVE THE DEVICE WHILE ADJUSTING.


Undo the height adjustment handle (see Fig. 4.1), set the tube height required...



...then tighten the height adjustment handle (see Fig. 4.2).

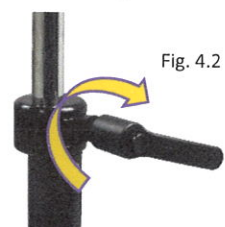

After tightening, the handle it can be repositioned to a more convenient position to minimise the risk of snagging, etc, without loosening the handle. Pull the handle outwards away from the handle body (Fig. 4.3), with the handle still pulled out turn to the required position (Fig. 4.4), then release the handle and allow it to re-engage with the handle body (Fig. 4.5).


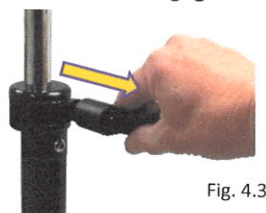


5. Adjusting the **SA-MF** Mini Floorstand Short Arm

The SA-MF features two Daessy L-Series Balljoints, an AFJ875 at the Post end and an AFJ750 at the Quick Release end. This procedure applies to both Balljoints. The tool required is the 5mm Allen Key supplied (and repositionable handle for the AFJ750 if preferred).

NEVER FORCE THE BALLJOINTS AND ALWAYS SUPPORT OR REMOVE THE DEVICE WHILE ADJUSTING

Loosen the bolt in the centre of the Balljoint (see Fig. 5.1). Set the position required and retighten the Balljoint (see Fig. 5.2). Do not overtighten, if using a Torque Wrench (such as our PTW12) set the torque to 12nm.

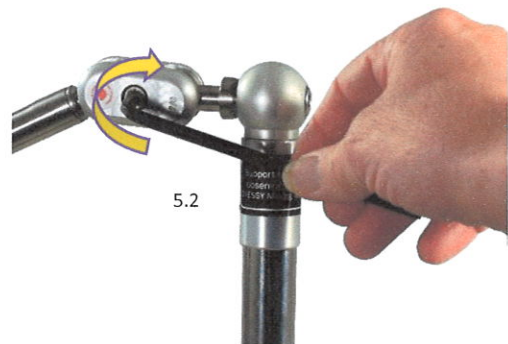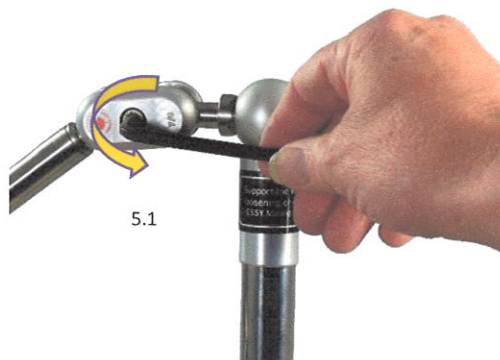

6. Adjusting the **PAQR** Mini Floorstand Post Articulating Quick Release Base

The PAQR is adjusted with a single Repositionable Handle.

ALWAYS SUPPORT OR REMOVE THE DEVICE WHILE ADJUSTING

Loosen the handle in the PAQR (see Fig. 6.1). Set the position required and retighten the handle (see Fig. 6.2).




Follow the same procedure as in section 4 to reposition the handle as required.

Looking after your Mini Floorstand

To clean the base and post, use warm water and a soft cloth or sponge.
Do not leave the base or post wet.



Do not use anything abrasive such as paper towels or scouring sponges to clean the base or post with.

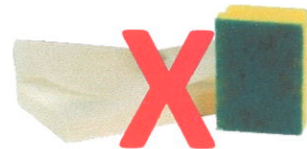

Polish as required but for the base and post we recommend regularly using a smooth and non-abrasive polish.



Do not use restoring polishes on the base or post, even if scratched.



If the post release collar becomes stiff or noisy when used, lubricate with Vaseline™ or similar. We recommend applying with a paint brush.



