



I-CLICK


Mounting example :
Spectacle frames

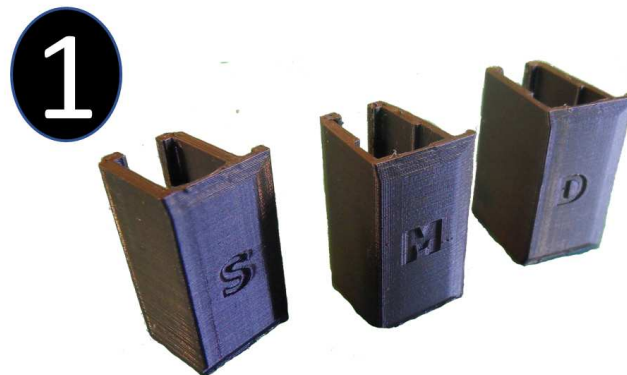

Choose the best size of *U Clamp* to suit the frame. Short, Medium or Deep

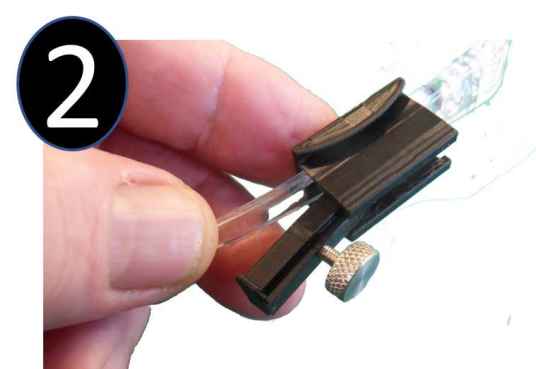


Wrap the *U Clamp* around frame and insert the thumb screw *Slide Assembly*. The wider T end should be at the end the rod will be from

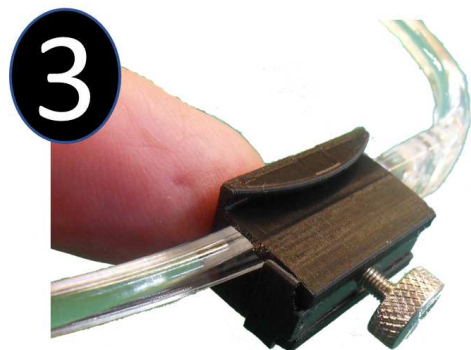


Fully insert the *Slide Assembly* with thumbscrew undone, leaving a gap for the support rod.


Important Tip
After you have found a good shape for the mounting rod that works, resist continuously bending the rod until you have tried repositioning whatever the rod terminates in.
ie In this example adjust the spectacles on the user

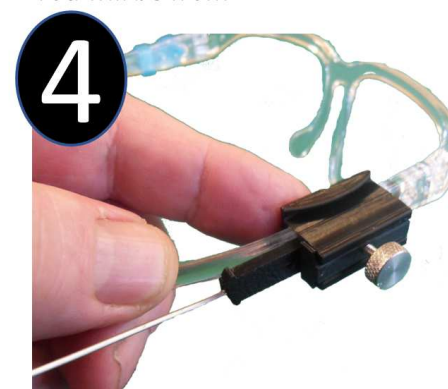

FULLY insert the magnetic *I-Click Mounting Rod (IMR)* between the frame leg and the Slide Assembly, taking care that the magnetic socket on the IMR for the sensor is facing the user.

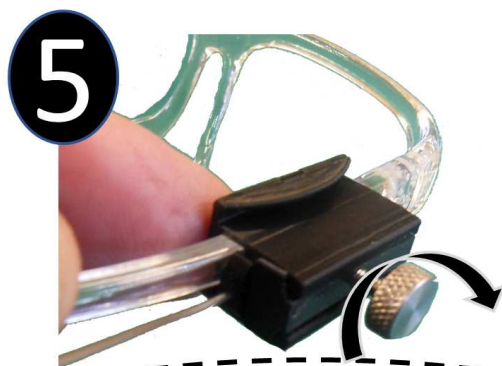

Tighten the thumbscrew to secure. **Do not overtighten the thumbscrew beyond that necessary to secure**

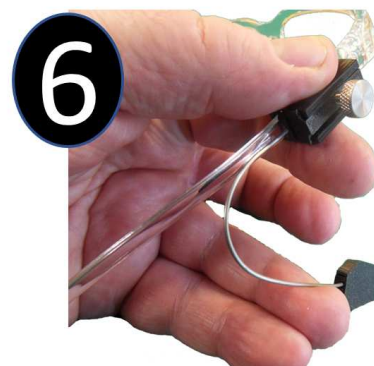

Without straining the leg of the spectacles bend and shape the wire into the approximate shape for use



To suit the user the IMR wire can be rear or front mounting and above or below the arm of the spectacles

Tip: For movement around the eye try the rear option first

