

BJ-254 Control USB M

Quick User Guide



www.bjliveat.com

1. Equipment description

Use **Control USB** to control any electronic device and elements of your environment (lights, blinds, doors, ceiling hoist, etc.) from your computer.

Control USB is comprised by the Control USB device and BJ Control software:


Software BJ Control



2



4

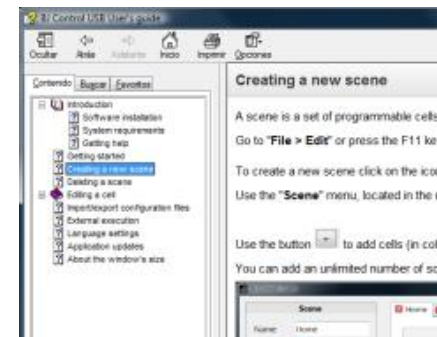

3

Customize the remote control interface and functionality with BJ Control software.

3. Proper care and use the device

- **Keep the device dry.** Moisture and liquids may damage the parts or electronic circuits of the device. Water damage to your device may void the warranty. **Dust may cause** the device to **malfunction**.
- **Do not drop the product** or subject it to impacts. It may be damaged if dropped. **Do not paint** your device or attach stickers to it.
- **Do not store the device near magnetic fields.**
- **Do not disassemble,** modify or repair the device. Any changes or modifications made to your device could void the warranty.
- When **cleaning** the device: **Use a towel or soft cloth** to clean the device and charger. **Do not use chemicals,** detergents or solvents.

Do not use it for any purpose other than that for which it was designed.

2. How to use

Download **BJ Control software** from:


<http://www.bjliveat.com>

1

Save the file in your computer and start the installation process.

The Setup Wizard will guide you through the installation process.

Plug the Control USB device into the USB port of your PC

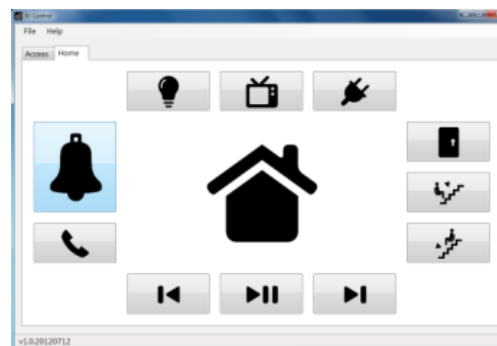

Check the **user's guide** within the **Help menu** to learn everything about this product.

Medical devices

Radio equipment can interfere with the functioning of medical devices that are not suitably protected. Consult a physician or the manufacturer of the medical device to determine whether it is adequately shielded from external radio frequency energy.

Implanted medical devices

Manufacturers of medical devices recommend a **minimum separation** of 15.3 centimetres (6 inches) between a wireless device and an implanted medical device such as a pacemaker or implanted cardioverter defibrillator, to avoid potential interference with the medical device.

Operating environment

- **Other products that use the same frequency band** (433 MHz) as this remote control **can interfere** with the control of the receivers. Make sure that neither the remote nor the receivers are located near these sources.

5. General information

Disposal

Waste electrical products should not be disposed of with household waste!

Dispose of the waste product via a collection point for electronic scrap or via your specialist dealer. Any unauthorized tampering with, or modifications to the product will annul this warranty.



Warranty

Within the statutory warranty period we undertake to rectify free of charge by repair or replacement any product defects arising from material or production faults.

Any unauthorized tampering with, or modifications to the product will annul this warranty.

Conformity

This product complies with the essential requirements of the Directive 2004/108/EC.

The Declaration of Conformity can be downloaded from the product's webpage at www.bjliveat.com.

Customer Service

If, despite of correct handling, faults or malfunctions occur or if the product was damaged, please contact us at the address below:

BJ Adaptaciones

C/ Mare de Déu del Coll, 70, Bajos
08023 Barcelona

España

Phone: +34 93 285 04 37

Fax: + 34 93 553 56 34

e-Mail: info@bjliveat.com

www.bjliveat.com

4. Technical features

- USB connection 2.0
- Power supply: 5V dc, 200mA
- Size: 50x25x13mm
- RF working frequency: 433MHz

BJ
Live!
smart solutions
that suit you